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 **TECNO-GAZ**
DENTAL AND MEDICAL EQUIPMENT

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TECNO-GAZ S.p.A.

TECNO-GAZ wishes to thank you for your purchasing one of its range of products and you are invited to read this booklet thoroughly.

This booklet contains all of the information necessary to use the product correctly.

Please adhere closely to the notices contained within and conserve the booklet, without making changes, in a place where it is easy to find and is accessible to the operator.

The contents of this manual can be changed, without warning or additional obligations, so that the manufacturer can include variations and improvements.

It is prohibited to copy or translate any of this booklet without the written permission of TECNO-GAZ.

Our Company remains at your complete disposal to provide any kind of information you may need for operating the appliance.

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*complies with the regulations in force regarding the Electromagnetic Compatibility Directive **89/336/EEC** and the Low Tension Directive **73/23/EEC** and their relative amendments and modifications.*

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POSITION: President of the Board of Directors

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For information regarding assistance services and retailers please contact your retailer or directly TECNO-GAZ Company:

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**P.S.: TECHNICAL ASSISTANCE IS DIRECTLY CARRIED OUT BY THE TECNO-GAZ
CONCESSIONAIRE.**

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WARRANTY

The device is covered by the warranty for a period of twelve months excluding electric components and those of the heating system which have a warranty period of six months. This warranty begins the date the device is delivered to the client, proved by the warranty coupon which has to be returned properly filled in, stamped and signed by the retailer. Repairing and replacement under warranty conditions of any components will be carried out exclusively by techno-gaz; the warranty does not cover travel expenses of personnel, neither packaging and transport expenses.

The warranty does not cover lamps, fuses, glassworks and the failures or damages caused by

wrong maintenance, negligence, unskilfulness or any other causes the constructor is not to be hold responsible for. Those materials subjected to regular wearing out are excluded from the warranty. TECNO-GAZ does not recognize the right to have the whole machine replaced.

The warranty does not meet the request of indemnity for direct or indirect damages, caused to goods and persons, provoked by the inefficiency of the device. TECNO-GAZ will not recognized any indemnity in case the machine is not used. The warranty will be automatically declined in case the machine is tampered, repaired or modified by the client or by third parties having not been authorized by tecno-gaz.

The purchaser is compelled to contact the retailer or after-sales service centres indicated by TECNO-GAZ for the repairing. The components replaced during the warranty period must be sent back to TECNO-GAZ porto franco. In case the replaced components are not be sent back to TECNO-GAZ the purchaser will be charged for them.

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INDEX OF PARAGRAPHS AND PICTURES

01	ECC CONFORMITY STATEMENT	3
02	EEC MARKING	4
03	LIST OF RETAILERS AND TECHNICAL ASSISTANCE SERVICES IN ITALY	4
04	WARRANTY	5
05	TECHNICAL FEATURES	21
06	INSTRUCTIONS FOR PACKAGING, STORAGE, TRANSPORT AND HOISTING	22
07	INTENDED USE OF EQUIPMENT	22
08	DESCRIPTION OF THE EQUIPMENT	23
09	EQUIPMENT'S PICTURE	24
10	INSTRUCTIONS FOR THE INSTALLATION	25
11	CONTROLS LIST	25
12	OPERATING INSTRUCTIONS	26
13	DANGEROUS SITUATIONS DURING OPERATIONS	28
14	HARMFUL SUBSTANCES EMITTED BY THE EQUIPEMENT	28
15	MANUAL OPERATIONS THAT MAY CAUSE DAMAGES	29
16	SAFETY EQUIPMENT FOR THE OPERATOR	29
17	ELECTRICAL SYSTEM DIAGRAM	30
18	LIST OF THE ELECTRICAL SYSTEM COMPONENTS	30
19	WARNING, SUGGESTIONS AND USEFUL INFORMATION	31
20	ISTRUCTIONS FOR THE EQUIPMENT DEMOLITION	32
21	WARNING, SUGGESTIONS AND USEFUL INFORMATION	32

PAR. 05**TECHNICAL FEATURES****MECHANICAL FEATURES**

HEIGHT	220mm
LENGHT	300mm
WIDTH	200mm
WEIGHT	6Kg
MATERIAL	STAINLESS STEEL AISI 304

USEFUL DIMENSIONS OF THE COMPARTMENT:

HEIGHT	75mm
LENGHT	220mm
WIDTH	120mm
WATER TRAY CAPACITY	1.6 Lt.

ELECTRICAL FEATURES

POWER SUPPLY	220V
FREQUENCY	50Hz
PHASES	1 + Neutral
FUSE CURRENT	4 A
ABSORBED POWER	0.7 KW
TYPE OF PROTECTION	PE
CLASS	1
SUPPLY CABLE mt.2	2x1.5mmq+T
CONNECTOR	16A+ ground

PAR. 06**INSTRUCTION FOR PACKAGING, STORAGE, TRANSPORT AND HOISTING**

The equipment DECERA does not need any special packaging since it is very small and its weight is very low.

The equipment is sent protected by a cellophane bag, wrapped in a foam polystyrene cover, placed within a corrugated cardboard box and sealed with adhesive tape.

ATTENTION: Do not abandon packaging material in the environment, recover all those materials that can be recycled and entrust any the other products to a company specialized in disposing of waste.

The packed equipment has to be kept in a dry place and at a temperature compressed between -10°C and +40°C.

It is allowed to pack the equipment one on top of other.

PAR. 07**INTENDED USE OF EQUIPMENT**

DECERA is a device used to melt the cylinders' was moulds; to be used only for dental technology.

It is not allowed to use the device for other purpose or in environment different from those the DECERA has been manufactured for.

ATTENTION.

The device DECERA does not melt thermoplastic stamps.

PAR. 08**DESCRIPTION OF THE EQUIPMENT**

The dewaxer is a very small device composed by two parts:

A) The piece of furniture, pos.51 pag.8.

B) The control board, pos. 52, pag. 8.

The piece of furniture houses a tray that contains the following equipment:

Pos. 56, thermal resistance's protection.

Pos. 57, wax-holder.

Pos. 58, cylinder-holders' grid.

Pos. 61, the cover.

A thermal resistance, pos.6 , located on the lower part of the tray , produces the necessary heat to boil the water.

The temperature is controlled by the bulb thermostat mounted over the resistance.

The control board contains the electrical equipment necessary to control the temperature (bulb thermostat) pos. 3-4-5, general switch pos. 1 and the emergency pilot light pos. 12.

PAR. 09

EQUIPMENT'S PICTURE

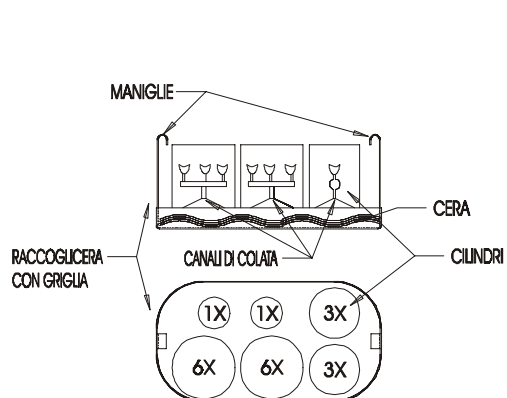
Fig.B


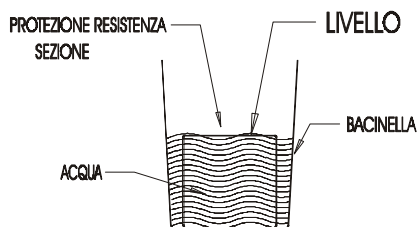
Fig.A

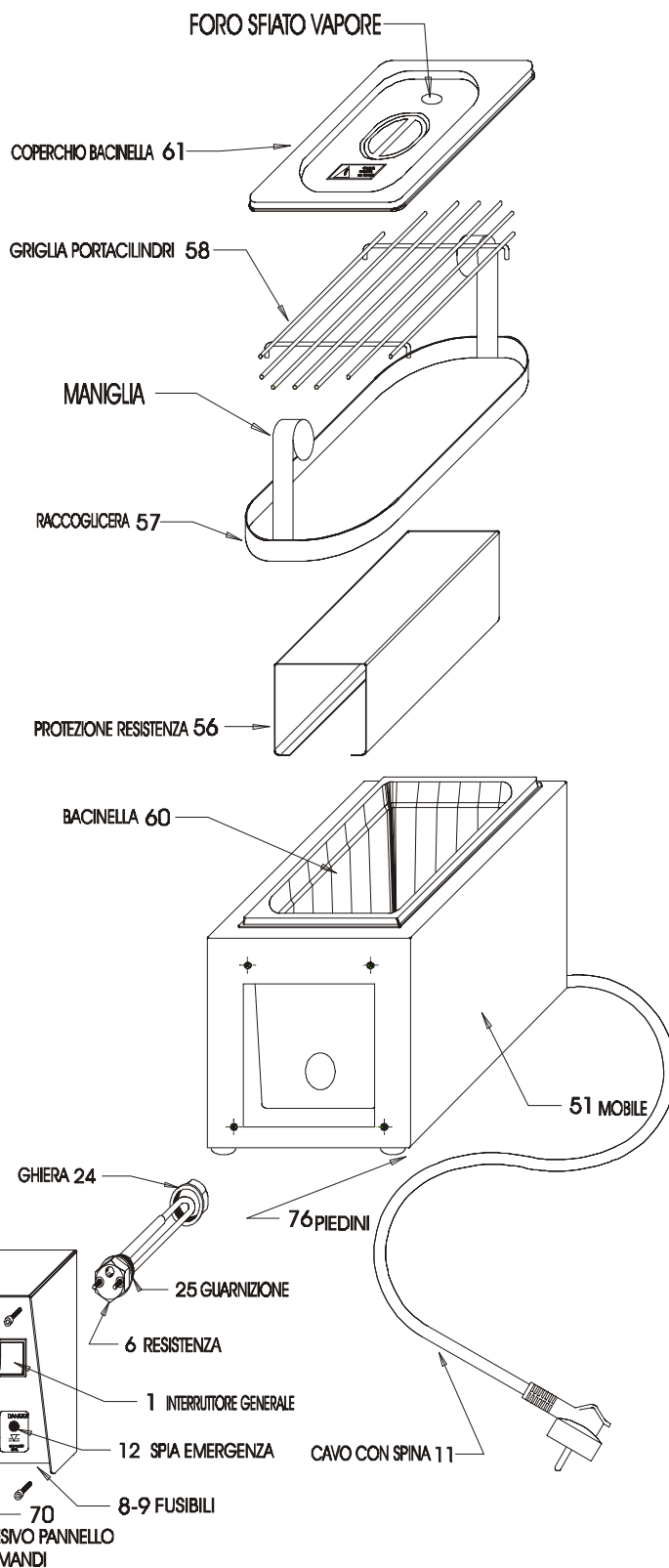PAR. 10**INSTRUCTION FOR THE INSTALLATION**

The dewaxer for cylinders "DECERA", has to be placed on an horizontal plane, under a suctioning hood or in a well ventilated area.

Connect the electric power plug pos. 11 to an alternate monophase outlet having a voltage of 220V 50Hz.

This outlet has to supply 1 KW power, it has also to be equipped with an earthing wire and a safety device.

PAR. 11**CONTROLS LIST**

THE CONTROL BOARD, Pos. 52, Pag. 8, INCLUDES THE FOLLOWING ELEMENTS:

- A) Pos.1 General switch.
- B) Pos.5 Thermostat for controlling and regulating the resistance's temperature.
- C) Pos.12 Pilot light that shows the thermostat is energized and therefore the presence of an emergency situation CAUSED BY THE LACK OF WATER INSIDE THE TRAY.

PAR. 12**OPERATING INSTRUCTIONS**

DECERA, as already described in par. 06 of pag. 6 is an equipment meant for melting cylinders' wax moulds.

To melt the wax it is necessary to heat the cylinders to such an steady temperature that allows the material to get melted without burning.

To obtain such a process, TECNO GAZ has adopted a steam heating system that allows the treatment to be carried out as described.

To dewax the cylinder it is necessary to carry out the following instructions:

A) Remove the cover pos. 61, pag. 8 and the wax-holder with the cylinder-holders grid pos. 57-58.

B) Fill in the tray, pos. 60, with tap water or demineralized water.

The level of the water admitted in the tray is indicated by the superior part of the resistance's protection, pos. 56; the tray has to be filled with water up to the above-mentioned surface. fig. A, pag. 8.

C) Arrange the cylinders on the grid pos. 58, lodged in the wax-holder pos. 57.

The cylinders have to be placed in vertical position, with the casting channels turned to the grid and should not come out from the wax-holder, fig. B, pag. 8.

D) Put the so composed group in the tray and close the cover; position the steam outlet hole to the opposite part of the control board.

E) Regulate, using the handle pos. 3, the thermostat to a temperature of 120°C.

The water have to boil continuously in order to produce the necessary steam needed to melt the wax.

F) Switch on the general switch, pos. 1, to start the treatment.

It takes DECERA more or less fifteen minutes to make the water boil.

The time of the treatment depends on the type of wax used, on the dimension of every single cylinder and specially on the diameter of the casting channels from which the wax flow out.

The fully loaded dewax time, with N° 2 cylinders 9X, is about 70 minutes, included the boiling time.

The fully loaded dewax time for cylinders 1X is reduced to 40 minutes, included the boiling time.

The above mentioned data is just indicative and does not represent any certainty applicable to the several pivoting methods.

Only the experience tells the operator which is the better treatment's time for every single case.

DECERA HAS a water autonomy of about 60 minutes, in optimal conditions.

If you need to extend the treatment time, it will be necessary to add water as described in par. 17 pag. 15, point H-I.

- G) Once the treatment has finished, switch off the general switch pos. 1 and wait some minutes for the steam to become condensation.
- H) Put on the gloves, remove the tray cover and take off the wax-holder using the special handles; see figure in pag. 8.
- I) Take away the cylinders from the grid, remove the wax that is deposited in the wax-holder and clean the equipment.

These operations have to be carried out at the end of every treatment.

- L) Control the level of water in the tray as described in point B); if necessary add water.
Replace the equipment in position and close the cover.
DECERA is ready to carry out another treatment.

For further information about the use and maintenance of the device, see par. 17 pag. 15.

PAR. 13**DANGEROUS SITUATIONS DURING THE FUNCTIONING****A) Interruption of the power supply.**

If power supply fails while the device is operating, DECERA goes off and get cold.

When the power supply is restored the equipment start working regularly, getting warm again.

B) Lack of water.

The lack of water or the lowering of the water level under the position of the resistance increases of the resistance's temperature.

The thermostat gets engaged, cuts off the power supply to the resistance and signals the danger situation by lighting the pilot light "DANGER" pos.2.

This safety situation is only temporary; no sooner the resistance's temperature decreases, falling under the value set on the thermostat, the resistance goes on again.

This situation will keep repeating till the correct level of water in the tray has been refilled.

If the above mentioned situation often happen, in time it could compromise the equipment's functioning and its working life.

The operator has continuously to control the device's functioning.

PAR. 14**HARMFUL SUBSTANCES EMITTED BY THE EQUIPMENT**

During operations steam and gas come out from the wax moulds contained in the cylinders.

As already described in the previous paragraphs, it is convenient to place and use DECERA under a suction hood.

TECNO-GAZ WILL NOT TO BE HOLD RESPONSABLE FOR DAMAGES TO PERSONS OR GOODS CAUSED BY NEGLIGENCE OR UNSKILFULNESS OF THE OPERATOR.

PAR. 15**MANUAL OPERATIONS THAT MAY CAUSE DAMAGES**

All the manual operations regarding the cylinders' dewaxing, can be the cause of accidents of thermal nature.

It is not possible to protect the operator from these risks as the use of protective barriers do not allow the use of the device.

To avoid or limit burnings, the operator must wear suitable protective equipment.

TECNO-GAZ WILL NOT TO BE HOLD RESPONSABLE FOR DAMAGES TO PERSONS OR GOODS DURING MANUAL OPERATIONS.

PAR. 16**SAFETY EQUIPMENT FOR THE OPERATOR**

The operator should wear the following safety equipment while using DECERA:


- a) Thermal gloves
- b) Thermal apron

These equipment is not supplied by TECNO-GAZ.

Please contact a specialized retail.

PAR. 17

ELECTRICAL SYSTEM DIAGRAM

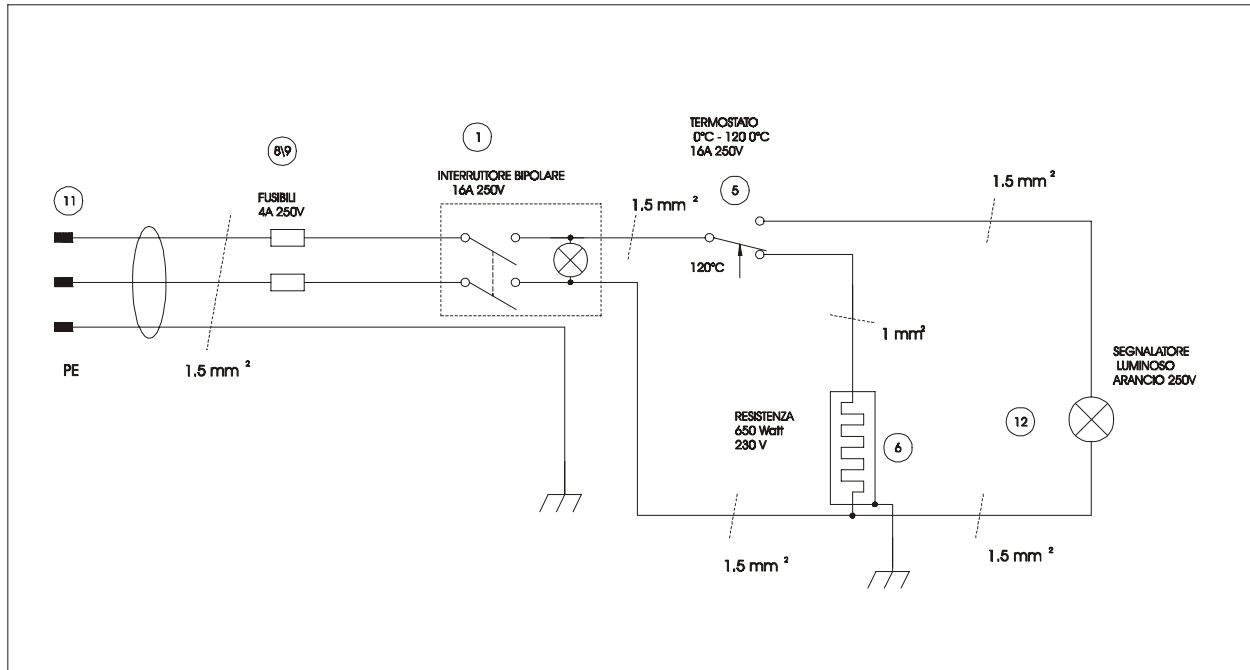PAR. 18

LIST OF THE ELECTRICAL SYSTEM COMPONENTS

Pos.	Q.ty	Description	Cod.
01	1	Bipolar platform balance switch, luminous, green, protected 250 V 10 A	CEQ0010
05	1	Bulb thermostat 16 A 250 V T 0°C / 120°C	CEB0002
06	1	Electric resistance 650 Watt 250 V	SDEA006
08	2	Board fuses	CEF0004
09	2	Fuses min. 5x20 4 A	CEF0009
10	1	Self blocking wire-inlet P16 ø 6 mm	CE02021
11	1	Cable with plug 3x1.5 mm²	CE06001
12	1	Orange luminous pilot light ø 10 mm 250 V	CEH0012

PAR. 19**WARNINGS, SUGGESTIONS AND USEFUL INFORMATION**

- A) Before starting operations, check the level of water and refill if necessary, see fig. a, pag. 8.
- B) Do not fill over the indicated level, since the water boils tumultuously, it will enter in the wax-holder and push the wax out.
- C) The total amount of water to fill in the tray is about 1.6 lt.
- D) Always close the cover.
- E) Avoid taking often off the cover, since this operation makes the cylinders cool down and the device uses more water so reducing its autonomy.
- F) Fill the tray only with water. Do not use liquids or solvents which getting heated up can provoke dangerous situations.
- G) During the treatment of the cylinders, the water contained in the tray, becomes steam and the level decreases. The time needed to pass from the maximum level to the minimum one (autonomy) is of about 60 minutes.
- H) When the sixty minutes of autonomy are over, and if you intend to continue with the treatment, it will be necessary to add immediately 1/2 lt. of water in the tray and repeat the refilling operation more or less every 40 minutes.
- I) To add water during operations, it is necessary to use a plastic container capable of containing the established quantity of water, and having a little hole (rubber ball). Lean the opening of the container on the wall of the tray, then press the container to make the water slide along the wall of the tray. Pay attention that the wax-holder does not touch the tray on the point where the water flows down.
- L) Leave always the thermostat regulated at the temperature of 120°C; if the thermostat is set to a lower temperature the water will not boil continuously, it will not produce steam and therefore the wax will not melt.
- M) The maximum height of the cylinders to be placed in the dewaxer is of 70 mm.
- N) The cylinders should not surpass the wax-holder limits as the melted wax will drop in the water and burn out when touching the resistance; specially when the level falls down.

PAR. 20**INSTRUCTIONS FOR THE EQUIPMENT DEMOLITION**

There are not particular instruction regarding the demolition and disposal of this equipment since stainless steel is quite recyclable and the remaining components have to be entrusted to the company specialized in disposing of these materials.

DO NOT LEAVE THE DEVICE IN UNATTENDED PLACES OR IN DANGEROUS POSITIONS.

PAR. 21**WARNINGS-SUGGESTIONS AND USEFUL INFORMATION**

- A) EVERYTIME YOU USE THE DEVICE CHECK THE WATER LEVEL, IF IS NOT ENOUGH REFILL IN THE WATER.
- B) PUT ON THE COVER.
- C) AVOID TAKING OFF OFTEN THE COVER, THIS OPERATION REDUCES THE AUTONOMY TIME OF THE DEVICE.
- D) THE MAX HEIGHT OF CILINDERS THE MACHINE CAN TREAT IN THE DEWAXER IS OF 60mm; CILINDERS OF 70-75mm CAN BEEN TREATED TAKING CARE TO PLACE THEM ON THE SIDE OF THE COVER'S HANDLE.

 **TECNO-GAZ** S.p.A.

DENTAL AND MEDICAL EQUIPMENT

Strada Cavalli, 4 - 43038 Sala Baganza, Parma, Italia

Tel. +39 0521 83.39.26 - Fax +39 0521 83.33.91

[http: www.tecnogaz.com](http://www.tecnogaz.com) - e-mail: info@tecnogaz.com

**COMPANY WITH QUALITY SYSTEM
CERTIFIED BY DNV**

== ISO 9002/EN 46002 ==